

Bioresorbable **PINS**

Bioabsorbable pin


The range of bioresorbable pins is mainly indicated for the fixation of small bone fractured in epiphyseal area.

- It is suitable for upper and lower bone extremities as well as pediatric indications.
- It is made of 100% poly(70/30;L/DL)lactide. This PLA follows controlled degradation process.
- Pins are stable during 8-10 months after implantation. Then pins are hydrolyzed in situ and totally replaced by bone after 24 months.

Bioresorbable PINS

Advantages


- Controlled degradation
- Synthetic & biocompatible
- Easy to use
- Mechanical strength
- Divisible

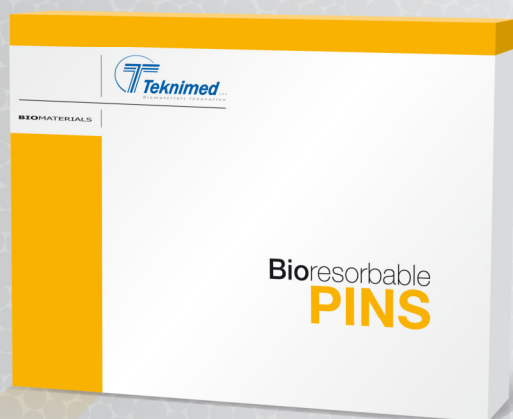

Resorbable pins are made of 100% poly (70/30;L/DL) lactide.

This PLA has a controlled kinetic degradation according to its initial molar mass. The upper the molar mass is, the lower the degradation is.

Our process guarantees that each absorbable pin has the same molar mass. Thus, all our pins present resorption at the same speed.

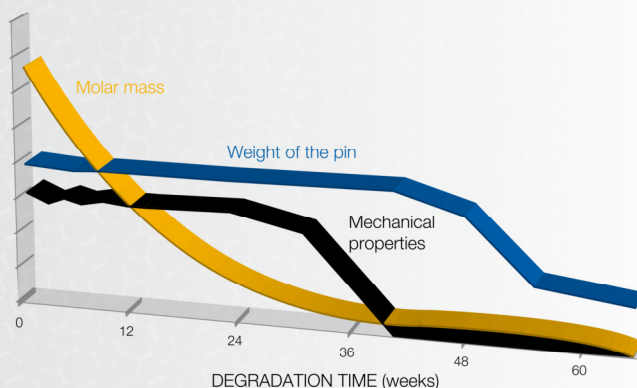

At 9 months, pins lose their mechanical properties by hydrolysis.

After 2 years, the pins are totally replaced by healthy bone.

Indications

Aligning and fixation of small bone fractured in epiphyseal area / Upper and lower extremities
/ Arthrodesis / Pediatric

Storage & products

Resorbable pins are delivered sterile and ready to use. The bioresorbable pin must be kept in its original unopened packaging in a clean and dry place at ambient temperature.

Each product is packed in double sealed bags and sterilised with ethylene oxide.

For any further information, please refer to the IFU.

DESIGNATION

BioresorbablePINS - Ø 2.0 mm x L 60 mm
BioresorbablePINS - Ø 2.4 mm x L 60 mm
 Kit **BioresorbablePINS** - Ø 2.0 mm x L 60 mm
 Kit **BioresorbablePINS** - Ø 2.4 mm x L 60 mm

REFERENCES

T7100220
 T7100224
 T7110220
 T7110224



Classe III



Head offices:

8, rue du Corps Franc-Pommiès
 65500 VIC en BIGORRE (France)
 Tél (33) 5 62 96 88 38
 Fax (33) 5 62 96 28 72

Administration and supply-chain

offices:
 Z.I de la Herray
 65500 VIC en BIGORRE (France)
 Tél (33) 5 62 96 88 38
 Fax (33) 5 62 96 28 72

Sales and engineering

departments:
 Z.I. de Montredon - 11, rue Apollo
 31240 L'UNION (France)
 Tél (33) 5 34 25 10 60
 Fax (33) 5 34 25 27 39

Manufactured by

